

NEW ASPIRET

DESK-TYPE

A surgical suction unit indicated for minor surgical applications and post-operative therapy.

Thanks to its compact size it finds application in out-patient and in-patient care, elderly care and in private care.

The device is suitable for nasal, oral and endotracheal aspiration of bodily fluids (mucus or catarrh) from adults and children.

MAIN FEATURES



Non-stop operation

Approved for non-stop operation, no interruption to cool down is required.



Small size - high performance

Powerful engine, reaching quickly a max vacuum of 0,75 Bar.



High-quality materials

The chassis is made of VO(UL94) ABS durable plastic, flame retardant, corrosion resistant and rust-proof.



Long lasting accessories

Made of resistant, autoclavable materials to guarantee safety and hygiene over time.



Easy to handle

Tiny and lightweight, can be placed anywhere and carried with ease.



Best Seller configuration

NEW ASPIRET +1000 mL reusable jar in Makrolon Polycarbonate (121°C)
Product code: RE 310001

Select your configuration at page 46

The warranty lapses immediately if the unit is used without all the required accessories or with non-original CA-MI accessories.



APPLICATIONS

HOME-CARE	ENT
ROOM / WARD	TRACHEOSTOMY

TECHNICAL DATA

MOTOR	Oilless and maintenance-free piston pump
POWER FEEDING	220-230 V / 50 Hz or 120-127 V / 60 Hz
ISO 10079-1 CLASSIFICATION	High vacuum / Low flow
MAX VACUUM adjustable	-0,75 Bar -75 kPa -563 mmHg value at sea level different altitudes may affect it
MAX FREE AIR FLOW RATE without accessories	15 L/min
NOISE LEVEL	59,6 dB(A)
DUTY CYCLE	Non-stop operation
IP PROTECTION DEGREE	IP21
WEIGHT	2,20 Kg
DIMENSIONS	19 x 23 x 15 cm
WARRANTY	2 years
SHIPPING CARTON	4 units
MANUFACTURING COUNTRY	Italy

CONVENIENT HANDLE



INCLUDED ACCESSORIES

REUSABLE JAR OR DISPOSABLE SYSTEM

depending on configuration

REUSABLE SILICONE TUBES

Ø 6 / 10 mm

BI-CONICAL CONNECTORS

Ø 8 / 9 / 10 mm

ANTIBACTERIAL & HYDROPHOBIC FILTER

CH20 SUCTION CANULA

POWER CORD WITH SCHUKO PLUG

other plugs available upon request

OPTIONAL ACCESSORIES

Reusable Jars in Makrolon®, APEC® Polycarbonate or Polysulfone

AVAILABLE SIZES:
1000 mL / 2000 mL

AVAILABLE SETS:

- › Complete set (Jar+cover with overflow system+ O-ring)
- › Jar (without cover)
- › Cover with overflow system
- › O-ring

Further details at page: 41



OPERATING WITH:

- REUSABLE COMPLETE TUBES' SET
- 3 silicone tubes (Ø 6/10mm, L: 140cm/ 13cm each)
- + 1 bi-conical connector Ø 8-9-10mm
- + 1 Antibacterial and hydrophobic filter


Disposable Flovac® System

AVAILABLE SIZES:
1000 mL / 2000 mL

AVAILABLE SETS:

- › Reusable Flovac® Canister
- › Disposable Flovac® Liner
- › Male Connectors for Flovac® Liner (5's pack)

Further details at page: 43

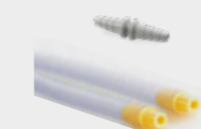

OPERATING WITH:

- FLOVAC COMPLETE TUBES' SET
- 2 tubes (Ø 6 / 10mm, L: 140cm / L: 26cm)
- + 1 bi-conical connector (Ø 8 / 9 / 10mm)

Rolls of silicone tube

AVAILABLE SIZES:

- › Ø 6 / 10mm, Length = 1m
- › Ø 6 / 10mm, Length = 10m
- › Ø 6 / 10mm, Length = 50m



Cannulas

AVAILABLE MODELS

- › Aspiration Probe CH20 (set of 10's)
- › Yankauer Handle Flat Tip With Hole
- › Yankauer Handle Crown Tip With Hole
- › Yankauer Tube L= 180cm



Stands

Stand on 5 Castor Ø 50mm with brakes



Carrying Bags

AVAILABLE SIZES:
1000 mL / 2000 mL



NEW 350 HOSPIVAC





ON CASTORS


A professional surgical suction unit indicated for any surgical application guaranteeing high performances and versatility.

Thanks to its state-of-the-art maintenance-free pump it finds application in any operating theatres and general surgery.

The device is suitable for nasal, oral and endotracheal aspiration of bodily fluids (mucus or catarrh), tissues or bones from adults and children.

MAIN FEATURES

-  **Non-stop operation**
In order to guarantee the best performance, engine can work without interruption 24/24h.
-  **Change-over system**
Allows operator to switch to the next jar without interrupting suction.
-  **Safety protection**
3 protection systems against internal reflux.
-  **Multi-Purpose Rail System**
Flexible rail that allows the rapid change of containers.

 **Best Seller configuration**
NEW HOSPIVAC 350 +2x 2000 mL reusable jar
in Makrolon Polycarbonate (121°C)
Product code: RE 410356

Select your configuration at page 46

The warranty lapses immediately if the unit is used without all the required accessories or with non-original CA-MI accessories.



APPLICATIONS

GENERAL SURGERY	NEUROSURGERY
GINECOLOGY	OBSTETRICS
OPERATING THEATRE	DENTAL PRACTICE

TECHNICAL DATA

MOTOR	Oilless and maintenance-free piston pump
POWER FEEDING	230 V / 50 Hz or 110-127 V / 60 Hz
ISO 10079-1 CLASSIFICATION	High vacuum / High flow
MAX VACUUM adjustable	-0,90 bar -90 kPa -675 mmHg value at sea level different altitudes may affect it
MAX FREE AIR FLOW RATE	60 L/min
NOISE LEVEL	51.7 dB(A)
DUTY CICLE	Non-stop operation
WEIGHT	13 Kg
DIMENSIONS	46 x 42 x 85 cm
WARRANTY	2 years
SHIPPING CARTON	1 unit
MANUFACTURING COUNTRY	Italy



INCLUDED ACCESSORIES

Basic Setup

2 REUSABLE JARS OR DISPOSABLE SYSTEM

depending on configuration

AUTOCLAVABLE SILICONE TUBES

ø 8/14 mm

ANTIBACTERIAL & HYDROPHOBIC FILTER

single-patient

SAFETY TRAP BOTTLE

220 mL

BI-CONICAL CONNECTOR

ø 10/11/12 mm

JAR ACCOMODATION RINGS

size depending on the model

AIR SUCTION INLET

POWER CORD WITH SCHUKO PLUG

FS Setup

ALL THE ABOVE-MENTIONED ACCESSORIES

+ FOOT-SWITCH CONTROL

with intermittent or continuous operation

Full Setup

ALL THE ABOVE-MENTIONED ACCESSORIES

+ CHANGE-OVER SYSTEM

from jar to jar by soft-touch keys

+ ANTIBACTERIAL & HYDROPHOBIC FILTER

single-patient

OPTIONAL ACCESSORIES

Reusable Jars in Makrolon®, APEC® Polycarbonate or Polysulfone

AVAILABLE SIZES:
2000 mL / 4000 mL / 5000 mL

AVAILABLE SETS:

- › Complete set (Jar+cover with overflow system+ O-ring)
- › Jar (without cover)
- › Cover with overflow system
- › O-ring

Further details at page: 41



OPERATING WITH:

REUSABLE COMPLETE TUBES' SET
3 silicone tubes (ø 8/14mm, L: 150cm/ 18cm each)
+ 1 bi-conical connector ø 10-11-12mm
+ 1 Antibacterial and hydrophobic filter




Disposable Flovac® System

AVAILABLE SIZES:
2000 mL / 3000 mL

AVAILABLE SETS:

- › Reusable Flovac® Canister
- › Disposable Flovac® Liner
- › Male Connectors for Flovac® Liner (5's pack)

Further details at page: 43

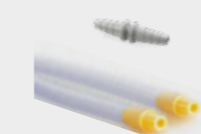

OPERATING WITH:

FLOVAC COMPLETE TUBES' SET
2 tubes (ø 8 / 14mm, L: 140cm / L: 26cm)
+ 1 bi-conical connector (ø 10/11/12mm)



Rolls of silicone tube

AVAILABLE SIZES:

- › ø 8 / 14mm, Length = 1m
- › ø 8 / 14mm, Length = 10m
- › ø 8 / 14mm, Length = 50m



Cannulas

AVAILABLE MODELS

- › Yankauer Handle Flat Tip With Hole
- › Yankauer Handle Crown Tip With Hole
- › Yankauer Tube L= 180cm

Foot-switch control

with plug-in system



Multi-purpose rail & accessories

Further details at page: 9



ASKIR C30

ON CASTORS

A surgical suction unit indicated for professional use.

Thanks to its powerful aspiration and to the double collection system it finds application in general surgery and

The device is suitable for nasal, oral and endotracheal aspiration of bodily fluids (mucus or catarrh) from adults and children.

MAIN FEATURES

Non-stop operation

In order to guarantee the best performance, engine can work without interruption 24/24h.



Maneuverable

Conveniently placed on a 5-castors stand to make its displacement easier.



High-quality materials

The chassis is made of UL94-V0 durable, flame retardant, corrosion resistant and rust-proof ABS plastic.



Long lasting accessories

All accessories are made in resistant, autoclavable materials to guarantee safety and hygiene over time.

Best Seller configuration

ASKIR C30 +2000 mL reusable jar in Makrolon Polycarbonate (121°C)
Product code: RE 410250

Select your configuration at page 46

The warranty lapses immediately if the unit is used without all the required accessories or with non-original CA-MI accessories.



APPLICATIONS

GENERAL SURGERY	NEUROSURGERY
GINECOLOGY	OBSTETRICS
OPERATING THEATRE	DENTAL PRACTICE

TECHNICAL DATA

MOTOR	Oiless and maintenance-free piston pump
POWER FEEDING	230 V / 50 Hz
ISO 10079-1 CLASSIFICATION	High vacuum / High flow
MAX VACCUM adjustable	-0,80 Bar -80 kPa -600 mmHg value at sea level different altitudes may affect it
FREE AIR FLOW RATE	40 L/min
NOISE LEVEL	61,5 dB(A)
DUTY CICLE	Non-stop operation
WEIGHT	6,20 Kg
DIMENSIONS	32 x 30 x 99 cm
WARRANTY	2 years
SHIPPING CARTON	1 unit
MANUFACTURING COUNTRY	Italy



INCLUDED ACCESSORIES

Basic Setup

2 REUSABLE JARS OR DISPOSABLE SYSTEM
depending on configuration

AUTOCLAVABLE SILICONE TUBES
ø 8/14 mm

ANTIBACTERIAL & HYDROPHOBIC FILTER
single-patient

BI-CONICAL CONNECTOR
ø 10/11/12 mm

POWER CORD WITH SCHUKO PLUG

FS Setup

ALL THE ABOVE-MENTIONED ACCESSORIES
+ FOOT-SWITCH CONTROL
with intermittent or continuous operation

OPTIONAL ACCESSORIES

Reusable Jars in Makrolon®, APEC® Polycarbonate or Polysulfone

AVAILABLE SIZES:
2000 mL

- AVAILABLE SETS:
- › Complete set (Jar+cover with overflow system+ O-ring)
 - › Jar (without cover)
 - › Cover with overflow system
 - › O-ring

Further details at page: 41




OPERATING WITH:
REUSABLE COMPLETE TUBES' SET
3 silicone tubes (Ø 8/14mm, L: 150cm/ 18cm each)
+ 1 bi-conical connector Ø 10-11-12mm
+ 1 Antibacterial and hydrophobic filter

Disposable Flovac® System

AVAILABLE SIZES:
2000 mL

- AVAILABLE SETS:
- › Reusable Flovac® Canister
 - › Disposable Flovac® Liner
 - › Male Connectors for Flovac® Liner (5's pack)

Further details at page: 43

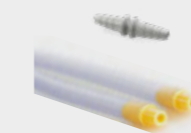

OPERATING WITH:
FLOVAC COMPLETE TUBES' SET
2 tubes (Ø 8 / 14mm, L: 140cm / L: 26cm)
+ 1 bi-conical connector (Ø 10/11/12mm)

Rolls of silicone tube

- AVAILABLE SIZES:
- › Ø 8 / 14mm, Length = 1m
 - › Ø 8 / 14mm, Length = 10m
 - › Ø 8 / 14mm, Length = 50m

Cannulas

AVAILABLE MODELS

- › Yankauer Handle Flat Tip With Hole
- › Yankauer Handle Crown Tip With Hole
- › Yankauer Tube L= 180cm

Foot-switch control

with plug-in system

