


# AKAI<sup>Med</sup>

## BABY INCUBATOR STANDARD CONFIGURATION

Main body(including the Transparent hood, Control system, Infant bed, Water trough), I.V.pole, Mattress, Drawers, Castors, Double wall hood.

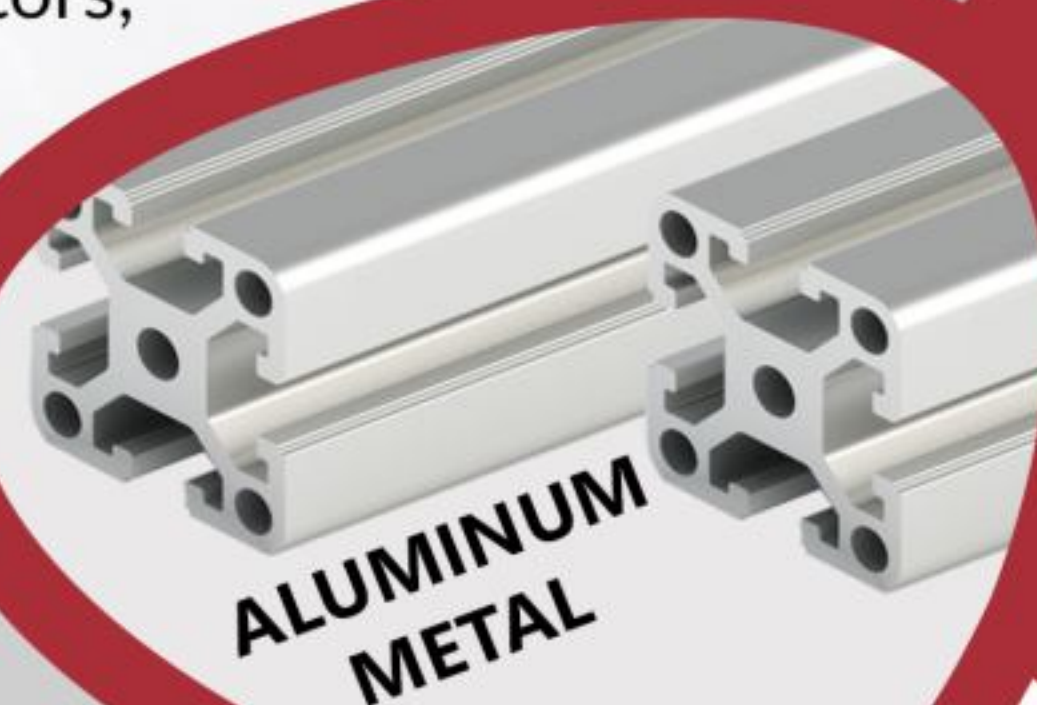

### PANDA 500 PLUS Plus Series





# Other features :

- Microprocessor based servo controlled temperature system.
- Control modes: Air mode and baby skin mode.
- Humidity is adjustable continuously.
- New Integrated controller for easy maintenance and more safety
- Set temperature, air temperature, baby skin temperature, timer and heating power are displayed separately by LED.
- Self test function and failure alarm.
- >37 °C temperature set function and could reach 38°C.
- Triple protection for over temperature with separate cut off device,more safety system.
- The inclination of infant bed is adjustable.
- Double wall hood, 4 operating windows and 2 iris ports.
- Independent locking device for front door.
- New design window cover, more easy operating
- With four drawers / IV pole and shelf.



## Technical data

Power supply	AC110/220V 60/50Hz
Power input	420VA

## Operating conditions:

Environment temperature	20°C ~ 30°C
Environment relative humidity	30%~75%RH
Environment air velocity of flow	<0.3m/s
Air Temperature control range	25°C~38°C
Skin temperature control range	32°C~38.5°C
Skin temperature sensor accuracy	±0.3°C
Temperature fluctuation	±0.5°C
Uniformity of mattress temperature	≤0.8°C
Warm-up time	≤45min

Internal noise level	≤55dB(A)
Infant bed tilt angle	±10°
Mattress size	65cm(L)*37cm(W)
Trough capacity	1200ml
Air filter	0.5um



## Standard configuration

Main body(including the Transparent hood, Control system, Infant bed, Water trough), I.V.pole, Mattress, Drawers, Castors, Double wall hood.


## Transport and storage

### Package

Each unit is packed in two wooden crate  
 \* (control unit and cabinet):112\*66\*99CM  
 N.W.:80KG G.W.:110KG  
 \* Hood: 98\*63\*75CM  
 N.W.:20KG G.W.:40KG

## Transport and storage

Environment Temp:-10°C--+70°C  
 Environment Relative Humidity: ≤80%  
 Atmospheric Pressure: 500~1060hPa

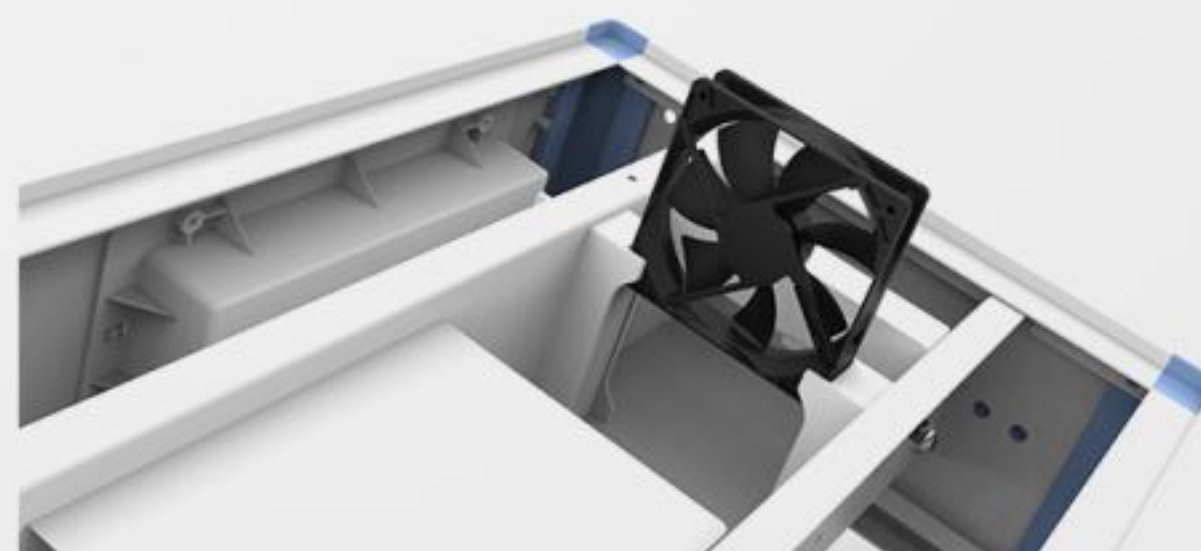

AKAIMED LIMITED CO.

Add: 706/F South seas centre,Tower 75, 2 Mody Road, Tsim Sha Tsui, Kowloon, Hong Kong

Email:international@akaimed.com

Tel:+852 9553 3913 /+852 5637 3743